


Quick start guide CCP3A-2217-EU

Care Call Pager



Any questions? Call us on free phone 0800 387 397

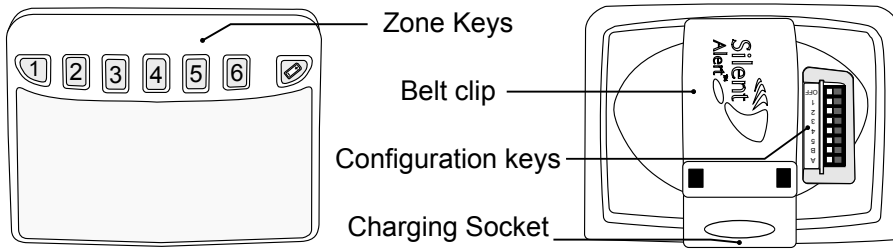

1.

Move the ON/OFF switch to the on position (to the left away from off).

The pager will illuminate each key, beep then vibrate. The pager is now ready to use.

2.

Install the trickle charger (Part no. CCT3A-2214-IN), insert the Pager into the charging slot and charge for 24 hours. The battery key on the Pager will light green briefly upon docking and the sounder will beep 3 times confirming the Pager is docked correctly.

If a pillow pad is fitted the charger will not beep and the pad will vibrate 3 times.

3.

Install the Care Call monitor units according to their instructions.

Remove the Pager from the Charger slot and test each monitor in turn.

We recommend a full 24 hour charge before using our equipment.

In an emergency a few hours charge may be sufficient for short time use but this is not advised.

4.

To re-set the Pager and stop the vibration, press and hold the key that is lit.

Each signal will automatically cancel after 20-40 seconds.

5.

The battery key will flash green every 45 seconds to indicate the pager is switched on.

The Pager must be recharged each night and if not in use for more than 24 hours should be switched off.

6.

If the Pager becomes flat it will vibrate and the battery key will light up. The battery key will flash orange until the Pager is recharged.

7.

If the low battery sequence alternates with another zone key the monitor configured to send that zone has a low battery. Replace the Monitor battery and transmit to test.

WARNING! Your device and its accessories are not toys. They may contain small parts. Keep them out of the reach of small children.

Any questions? Call us on free phone 0800 387 397

Increasing system range

It is possible to make the Pager more sensitive when larger areas need to be covered.

1.

With the Pager switched on, move switch 5 to the left (away from the belt clip).

2.

Operating Range

This will increase the range of the Pager from 100M (in open air), to approximately 1000M (in open air).

Vibration Alert

Users can choose between the standard audio alert or a discreet vibration pattern for each event.

1.

With the Pager switched on, move switch A to the left (away from the belt clip).

2.

The Pager will now vibrate when a signal is received.

This can be especially useful in the evening, effectively putting the Pager into "silent" mode.

Changing the system channel code

NOTE. In most cases it is not necessary to change the system code. However, when one or more systems are in close proximity, system codes can be used to avoid interference from other Care Call systems using up to a maximum of 16 channels.


The Pager is supplied with no system code switches set (system code 1). For reference this is the factory setting should you need to re-set the unit.

System codes can be set using key switches 1 - 4.

The diagram to the right shows the 16 possible combinations.

Be sure that the same system code is set on all monitors to be used in the system.

NOTE. Care Call receivers have an ON/OFF switch, Care Call monitors do not.

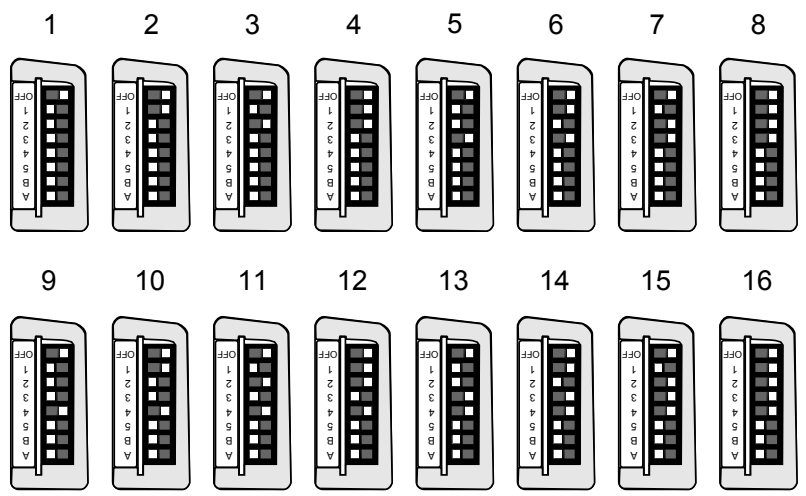

CORRECT

Here we have a Care Call Pager and Mini Monitor set to system channel 2.

Ensure when setting system channel keys you remember to ignore the Pagers **ON/OFF** switch or the equipment will be on different system channels .

The diagram to the right shows an incorrect attempt to set the system to channel code 2 (shown as an example to avoid possible confusion). The correct settings are shown to the left.

INCORRECT